


CONCIERGO Guiding
CATÉTER GUÍA Excellence™

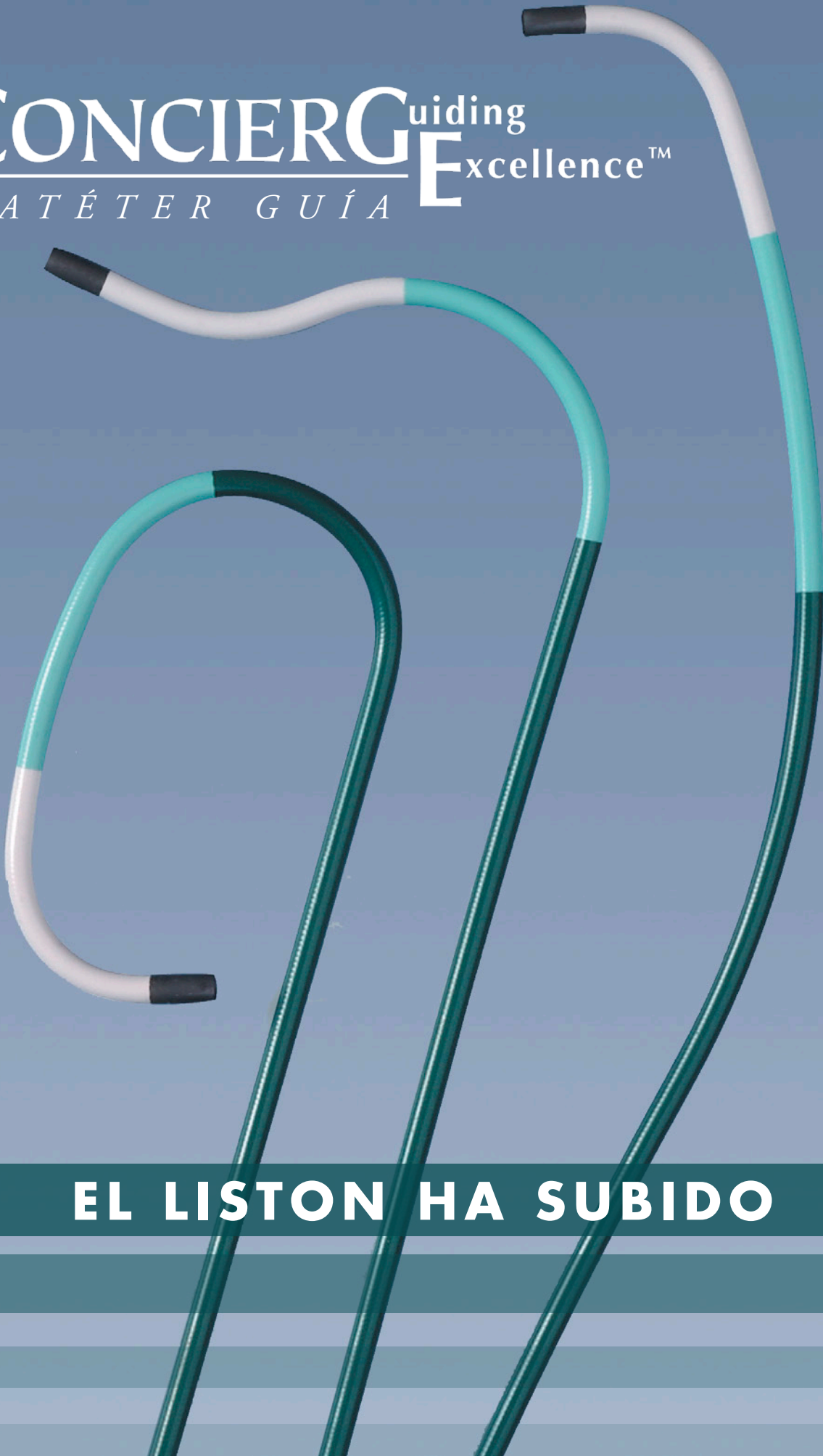

EL LISTON HA SUBIDO


CONCIERGE^{Guiding Excellence™}

En un estudio independiente, aproximadamente el **40%** de los cardiólogos intervencionistas afirmaron que su catéter de guía actual no satisface sus necesidades de **resistencia al acodamiento**.

Más del **80%** afirmaron que un **buen apoyo** o una **punta no traumática** son las características más importantes en un catéter guía.¹


El ConcierGE es el primer catéter que ofrece que ofrece la excelencia como catéter guía y que proporciona un rendimiento superior en las características que los hemodinamistas más valoran: resistencia al acodamiento, capacidad de apoyo y torque 1:1.

Resistencia a al acodamiento

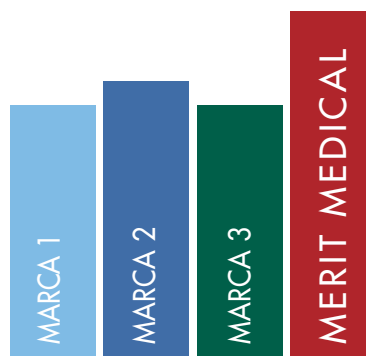

El ConcierGE tiene aproximadamente un **30% más de resistencia al acodamiento*** que las marcas líderes, lo que facilita el avance y la retirada de los dispositivos.

Capacidad de apoyo

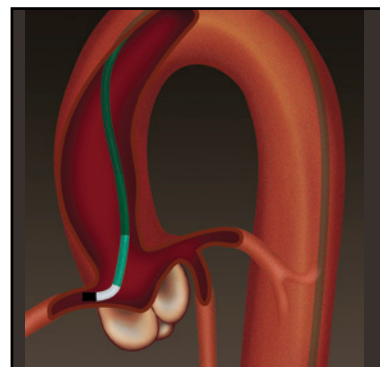

El catéter ConcierGE proporciona un apoyo fiable gracias **a la mezcla de Nylon utilizada en la construcción de su cuerpo y a la tecnología utilizada** en la malla interna del catéter. Esta característica permite avanzar dispositivos a la zona deseada sin comprometer la resistencia del cuerpo del catéter.

DESCRIPCIÓN DEL PRODUCTO			
5F	6F	7F	8F
DIÁMETRO EXTERNO DE DEL CATÉTER GUÍA			
0,070"(1,7 mm)	0,082"(2,0 mm)	0,092"(2,3 mm)	0,105"(2,7 mm)
DIÁMETRO INTERNO DE DEL CATÉTER GUÍA			
0,057"(1,4 mm)	0,070"(1,8 mm)	0,078"(2,0 mm)	0,088"(2,2 mm)
Malla - Acero inoxidable			
Capa interna - PTFE (politetrafluoroetileno)			
Capa externa - (mezcla propia de polímeros)			


¹ HRA: HealthCare Research & Analytics. Datos en archivo.

* Los resultados de las pruebas están basados en la comparación de un catéter guía de 6F. Datos en archivo.

** A 1/2 giro del conector del catéter. Datos en archivo.


EL LISTON HA SUBIDO

Punta no atraumática

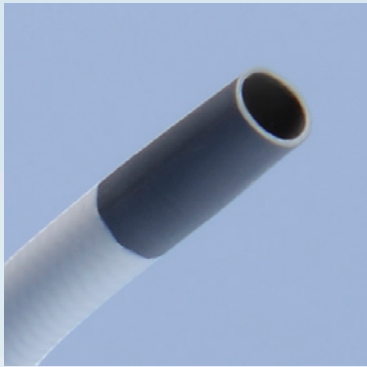


El ConcierGE tiene una de las **puntas más blandas**, lo que proporciona una navegación e intubación seguras de manera que médicos con distinto nivel de experiencia puedan trabajar con seguridad en anatomías complejas.

Torque 1:1

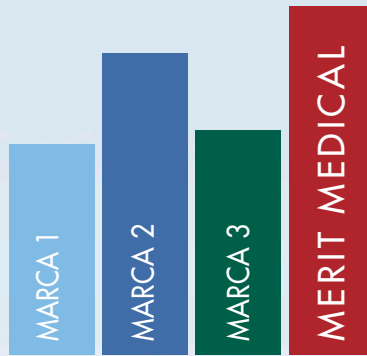


Los resultados de las pruebas del ConcierGE demuestran que tiene **un mejor torque 1:1**** que las marcas líderes, lo que ayuda a su intubación sin movimientos bruscos del catéter. *

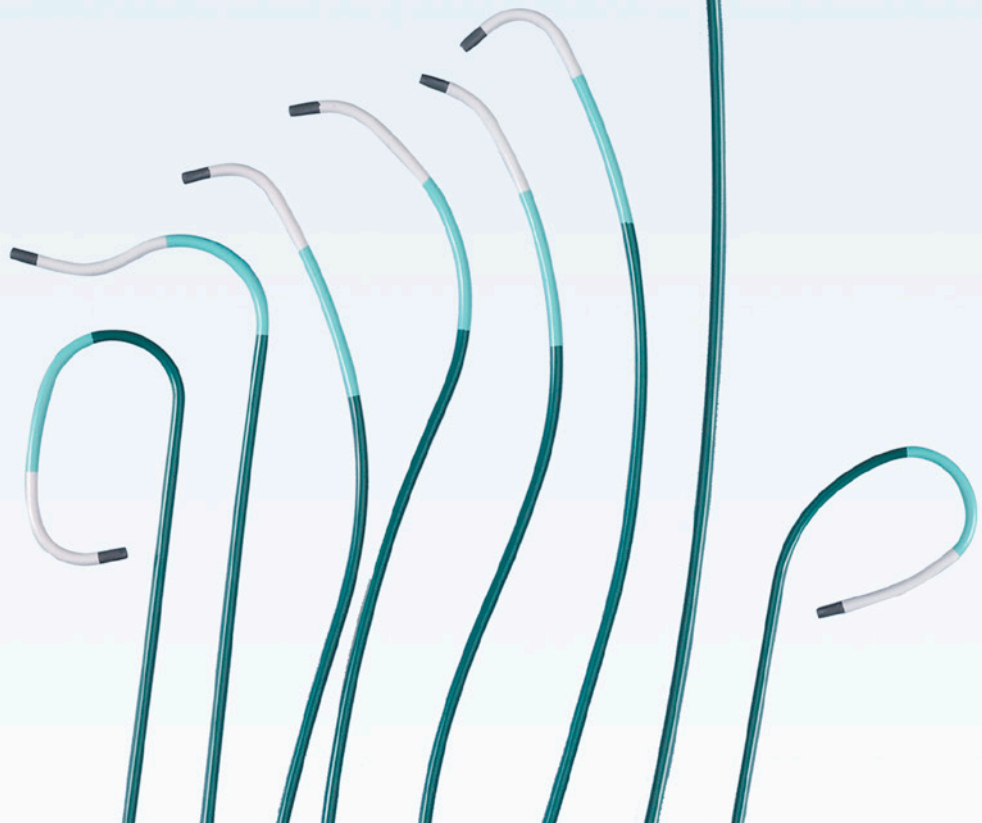


Números de catálogo de ConcierGE según su tamaño en French

Curva	Agujeros laterales	Punta	5F	6F	7F	8F
Judkins derecho	-	3,0	CGC5100JR30	CGC6100JR30	CGC7100JR30	CGC8100JR30
Judkins derecho	-	3,5	CGC5100JR35	CGC6100JR35	CGC7100JR35	CGC8100JR35
Judkins derecho	2	3,5	CGC5100JR35SH	CGC6100JR35SH	CGC7100JR35SH	CGC8100JR35SH
Judkins derecho	-	4,0	CGC5100JR40	CGC6100JR40	CGC7100JR40	CGC8100JR40
Judkins derecho	2	4,0	CGC5100JR40SH	CGC6100JR40SH	CGC7100JR40SH	CGC8100JR40SH
Judkins derecho	-	5,0	CGC5100JR50	CGC6100JR50	CGC7100JR50	CGC8100JR50
Judkins derecho	2	5,0	CGC5100JR50SH	CGC6100JR50SH	CGC7100JR50SH	CGC8100JR50SH
Judkins derecho	-	6,0	CGC5100JR60	CGC6100JR60	CGC7100JR60	CGC8100JR60
Judkins izquierdo	-	3,0	CGC5100JL30	CGC6100JL30	CGC7100JL30	CGC8100JL30
Judkins izquierdo	-	3,5	CGC5100JL35	CGC6100JL35	CGC7100JL35	CGC8100JL35
Judkins izquierdo	2	3,5	CGC5100JL35SH	CGC6100JL35SH	CGC7100JL35SH	CGC8100JL35SH
Judkins izquierdo (ST)	-	3,5	CGC5100JL35ST	CGC6100JL35ST	CGC7100JL35ST	CGC8100JL35ST
Judkins izquierdo	-	4,0	CGC5100JL40	CGC6100JL40	CGC7100JL40	CGC8100JL40
Judkins izquierdo	2	4,0	CGC5100JL40SH	CGC6100JL40SH	CGC7100JL40SH	CGC8100JL40SH
Judkins izquierdo (ST)	-	4,0	CGC5100JL40ST	CGC6100JL40ST	CGC7100JL40ST	CGC8100JL40ST
Judkins izquierdo	-	4,5	CGC5100JL45	CGC6100JL45	CGC7100JL45	CGC8100JL45
Judkins izquierdo	2	4,5	CGC5100JL45SH	CGC6100JL45SH	CGC7100JL45SH	CGC8100JL45SH
Judkins izquierdo	-	5,0	CGC5100JL50	CGC6100JL50	CGC7100JL50	CGC8100JL50
Judkins izquierdo	2	5,0	CGC5100JL50SH	CGC6100JL50SH	CGC7100JL50SH	CGC8100JL50SH
SBS 3,0	-	3,0	CGC5100SBS30	CGC6100SBS30	CGC7100SBS30	CGC8100SBS30
SBS 3,0	2	3,0	CGC5100SBS30SH	CGC6100SBS30SH	CGC7100SBS30SH	CGC8100SBS30SH
SBS 3,5	-	3,5	CGC5100SBS35	CGC6100SBS35	CGC7100SBS35	CGC8100SBS35
SBS 3,5	2	3,5	CGC5100SBS35SH	CGC6100SBS35SH	CGC7100SBS35SH	CGC8100SBS35SH
SBS 3,5 (Mod)	-	3,5	CGC5100SBS35-M	CGC6100SBS35-M	CGC7100SBS35-M	CGC8100SBS35-M
SBS 4,0	-	4,0	CGC5100SBS40	CGC6100SBS40	CGC7100SBS40	CGC8100SBS40
SBS 4,0	2	4,0	CGC5100SBS40SH	CGC6100SBS40SH	CGC7100SBS40SH	CGC8100SBS40SH
SBS 4,0 (Mod)	-	4,0	CGC5100SBS40-M	CGC6100SBS40-M	CGC7100SBS40-M	CGC8100SBS40-M

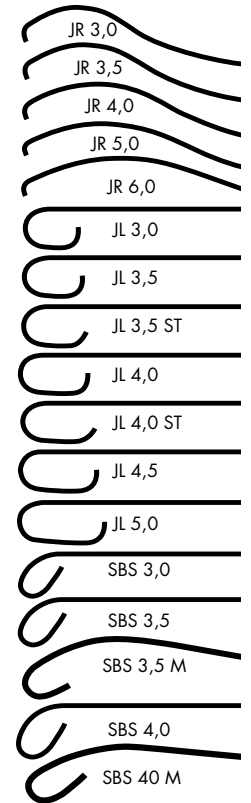

Números de catálogo de ConcierGE según su tamaño en French

Curva	Agujeros laterales	Punta	5F	6F	7F	8F
Amplatz derecho 1	-	AR1	CGC5100AR1	CGC6100AR1	CGC7100AR1	CGC8100AR1
Amplatz derecho 1	2	AR1	CGC5100AR1SH	CGC6100AR1SH	CGC7100AR1SH	CGC8100AR1SH
Amplatz derecho 2	-	AR2	CGC5100AR2	CGC6100AR2	CGC7100AR2	CGC8100AR2
Amplatz derecho 2	2	AR2	CGC5100AR2SH	CGC6100AR2SH	CGC7100AR2SH	CGC8100AR2SH
Amplatz derecho 3	-	AR3	CGC5100AR3	CGC6100AR3	CGC7100AR3	CGC8100AR3
Amplatz derecho 1	-	AL1	CGC5100AL1	CGC6100AL1	CGC7100AL1	CGC8100AL1
Amplatz derecho 1	2	AL1	CGC5100AL1SH	CGC6100AL1SH	CGC7100AL1SH	CGC8100AL1SH
Amplatz derecho 2	-	AL2	CGC5100AL2	CGC6100AL2	CGC7100AL2	CGC8100AL2
Amplatz derecho 2	2	AL2	CGC5100AL2SH	CGC6100AL2SH	CGC7100AL2SH	CGC8100AL2SH
Amplatz derecho 3	-	AL3	CGC5100AL3	CGC6100AL3	CGC7100AL3	CGC8100AL3
Multipropósito 1	-	MP	CGC5100MP1	CGC6100MP1	CGC7100MP1	CGC8100MP1
Multipropósito 1	2	MP	CGC5100MP1SH	CGC6100MP1SH	CGC7100MP1SH	CGC8100MP1SH
Bypass coronario izquierdo	-	LCB	CGC5100LCB	CGC6100LCB	CGC7100LCB	CGC8100LCB
Bypass coronario izquierdo	2	LCB	CGC5100LCBBSH	CGC6100LCBBSH	CGC7100LCBBSH	CGC8100LCBBSH
Bypass coronario derecho	-	RCB	CGC5100RCB	CGC6100RCB	CGC7100RCB	CGC8100RCB
Bypass coronario derecho	2	RCB	CGC5100RCBBSH	CGC6100RCBBSH	CGC7100RCBBSH	CGC8100RCBBSH
Mamaria interna	-	IM	CGC5100IM	CGC6100IM	CGC7100IM	CGC8100IM
Mamaria interna	2	IM	CGC5100IMSH	CGC6100IMSH	CGC7100IMSH	CGC8100IMSH
El Gamal	-	EG	CGC5100EG	CGC6100EG	CGC7100EG	CGC8100EG
Hockey Stick	-	HS	CGC5100HS	CGC6100HS	CGC7100HS	CGC8100HS
Ultimate 1	-	ULT1	CGC5100ULT1	CGC6100ULT1	CGC7100ULT1	-
Ultimate 2	-	ULT2	CGC5100ULT2	CGC6100ULT2	CGC7100ULT2	-

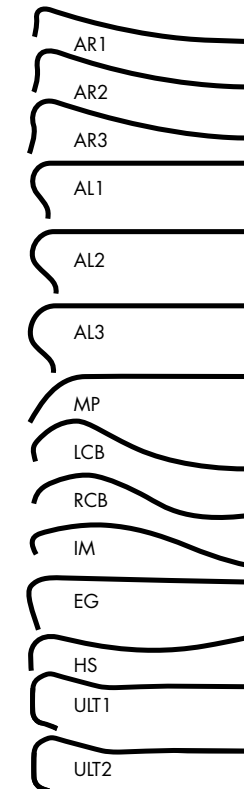

Todos los catéteres están disponibles con longitud de 100 cm

